

PRODUCT REGISTRATION

For your convenience you may register your Sonic Spatula™ one of three ways:

1. Visit www.TEISpa.com and click on "Register Your Product" under the Customer Care tab and submit your information
2. Email the information below to info@TEISpa.com
3. Mail the information below to TEI Spa. (Please visit www.TEISpa.com and click on "Contact Us" for an updated address)

Name _____

Company _____

Address _____

City _____

State _____ Zip _____

Telephone _____

Email _____

Date Purchased _____

Serial # _____

Store or Sales Rep _____

Limited 1 Year Warranty

If the Sonic Spatula™ tool fails to function within the specified warranty period because of defects in the materials or workmanship of the product, the tool should be returned to TEI Spa to be either repaired or replaced. This warranty does not cover any product that has been damaged by dropping, tampering, inappropriate usage, water damage or use with a power outlet other than the one provided. Proof of purchase required.

Experience the Benefits of
a Professional Spa Facial in
the Comfort of Your Home

A NEW VIBE IN CLEANSING

Push out impurities and clear away excess oil with the sonic vibration of the Sonic Spatula™. As the skin vibrates, the tissue and pores become softened, extracting environmental toxins, oil and makeup. On the flip side, the spatula aids in the absorption of water or gel-based skin products with the same sonic push, giving you double action in one tool. The energy of the ultrasonic waves brightens and evens skin tone while providing long-lasting results and radiant skin, all in just minutes a day.

- Vibrating ultrasonic waves exfoliate dead skin and extract dirt and oil from clogged pores
- Helps brighten and even skin tone
- Spatula shaped tip can be used to penetrate serums and gels for better efficacy
- Cordless & rechargeable
- Dual function in one tool



Phone: 310.231.0419 Fax: 310.231.0420

Sonic Spatula™: Designed by TEI Spa in California. Made in PRC.
HydraSkin Mask: Designed by TEI Spa in California. Made in Taiwan.
Lavender Hydra Mist: Designed by TEI Spa in California. Manufactured in USA.
© 2011 Tip Essentials, Inc.

TEI SPA ESSENTIALS

Sonic Spatula™

Ultrasonic Cleansing Tool



5% of all proceeds benefit the
Downtown Women's Center of Los Angeles



PRODUCT INSTRUCTIONS & SPECIFICATIONS

PATENT APPROVED 2011

Recommended & Used by Professional Estheticians

www.TEISpa.com

For All Skin Types

DIRECTIONS:

For Cleansing Mode:

1. Remove Protective Cap
2. On clean skin, spray included Lavender Hydra Mist onto treatment area to prepare for exfoliation
3. Press and hold the ON button for 3-5 seconds to turn device on. Red light indicates that Cleansing Mode is on.
4. With the spatula facing downward at an angle (as shown in Push-Out image below) gently exfoliate the pores of dead skin and impurities. Note: Treatment area must be moistened before use.
5. With a small amount of pressure on the skin, use the spatula in an upward and outward motion in small strokes.
6. Press and hold the ON button for 3-5 seconds to turn device off.

For Penetration Mode:

1. Remove Protective Cap
2. Apply preferred serum or gel on to treatment area to prepare for penetration.
3. Press and hold the ON button for 3-5 seconds to turn to cleansing mode. Press and hold the ON button **again** to turn to penetration mode. Green light indicates Penetration Mode is on.
4. With the spatula facing upward (as shown in Push-In image below) gently press the device flat on the skin and rub products into the skin in small circular motions.
5. Press and hold the ON button for 3-5 seconds to turn device off.


Note: Use the Sonic Spatula™ in Penetration Mode to deeply infuse the Hydra Skin's bio-cellulose content into the pores. See Hydra Skin Table for additional information.

For Sensitive Skin: Intermittent button can be used on both Cleansing Mode and Penetration Mode to lessen the intensity. Light will blink when intermittent mode is in use.



Push-Out impurities & exfoliate dead skin.

● Red light indicates CLEANSING MODE is ON

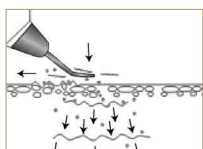

Push-In active ingredients in gels and serums.

● Green light indicates PENETRATION MODE is ON

RECOMMENDED USE:

Skin Type	Use Duration	Frequency
Oily Skin	3-4 minutes	3-4 times per week
Combination Skin	2-3 minutes	2-3 times per week
Dry/Sensitive Skin	1-2 minutes	1-2 times per week

NOTE: Daily use of Penetration Mode for enhancement of serums and moisturizers is recommended for all skin types.

HYDRA SKIN BIO-CELLULOSE FACIAL MASK

TEI Spa's Hydra Skin Facial Mask is made of bio cellulose, designed to increase moisture content, smooth fine lines and reduce pore size.

Directions for Use:

1. Wash face with gentle cleanser.
2. Carefully remove contents and unfold.
3. Remove outer layers to expose the Hydra Skin Facial Mask. **Note:** the mask is the middle layer, between the plastic and fabric layers.
4. Apply mask onto the face and adjust the position for comfort. **Note:** Use fabric layer on neck and chest for added moisture.
5. For optimum absorption of mask use the Sonic Spatula on Penetration Mode to deeply imbed Hydra Skin's bio-cellulose content into the pores.
6. Relax with mask on face for 15-20 minutes to allow for optimum absorption.
7. Remove mask when T-Zone has dried or after allotted time, you may apply water for easy removal.
8. Do not rinse face after removal of mask. Leaving mask residue on the skin allows for the fiber structure of Hydra Skin to transfer deeply into the pores.
9. Follow with daily skincare procedures.



SPECIFICATIONS:

Important Information about your unit:

When light above the ON button is blinking, the unit is set on Intermittent mode.

When recharge light on base is green, your skin spatula is fully charged and ready for use. Unplug base unit after use.

When recharge light on base is red, your skin spatula is being charged.

The Sonic Spatula™ Device will automatically power off after 10 minutes.

Clean the tool after each use with soap and a warm towel. If using the unit on multiple individuals, sterilize the spatula tip with alcohol after each use.

CAUTION:

Do not use if you:

- Have a pacemaker or other electrical medical devices
- Are pregnant
- Have recently had a cosmetic procedure
- Have recently had an acid or chemical peel or extreme sunburn
- Have skin cancer or severe acne

KEEP OUT OF REACH OF CHILDREN. DO NOT SUBMERGE IN WATER



When green light above the intermittent button is on, unit must be recharged.

THE SONIC SPATULA™ SET CONTAINS:

Sonic Spatula™
Protective Over-Cap
Lavender Hydra Mist
Hydra Skin Facial Mask
Gold Clutch
Recharge Base Unit
Power Adapter
Instruction Manual

Input recharge voltage DC12V
Input recharge power 3.6W
Frequency 15 - 35 KHz per second
Net weight Model: 112g Base: 98g
Size (mm) (L) 51.2 X (W) 26X (H)192 (R) 78.7X (H) 66
Stand by time: 720 Hrs
Charge time: 4 - 4.5 Hrs

www.TEISpa.com